

PUR Tracheostomy Tubes



FreeVent long term/permanent tracheostomy tubes are designed to provide the user with a safe, comfortable and flexible tube that is easy to handle.

FreeVent tubes come in a wide variety of sizes, with and without fenestration, speaking valves etc. The tubes are light, soft, pliable and harmless to tissue and adapt easily to the trachea.


The conical shape makes it much easier to insert and center the distal tip in the middle of trachea.

FreeVent tubes are made from tested medical grade Polyurethane that is thermosensitive i.e. softens in body temperature which increases the comfort for the patient.

Being transparent makes the FreeVent tubes more discreet and you can see through the cannula.

FreeVent PUR tubes do not need to be removed for X-ray or cobalt therapy, only the inner tube with silver speaking valve needs to be detached.

A low profile fitting on the FreeVent tube enables the user to comfortably wear tracheostomy HMEs and protectors.

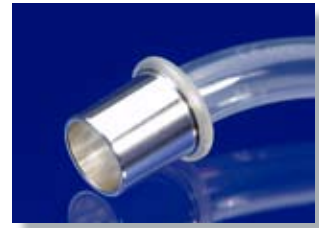

Plain inner cannula

FreeVent PUR Basic Tracheostomy Set consists of:

- 1 outer tube
- 1 plain inner cannula
- 1 inner cannula with 15 mm ISO connector

The PUR tracheostomy tubes provide a direct tracheal access for long term airway management. Thermosensitive PUR – provides sufficient rigidity for insertion but softens at body temperature which increases the comfort for the patient. The unique conical shape also makes it much easier to insert.

Feature	Benefit
Thermoplastic material	Soft and flexible at body temperature
Conical shape	Easier insertion & comfort
Wide range of models & sizes	Fits all
Low weight	High comfort
Excellent fit between inner and outer tube	No leakage
Low profile neck plate	Discrete design

Order information:

Product name	Size	Outer tube		Inner cannula		Length in mm	REF
		OD in mm	ID in mm	OD in mm	ID in mm		
FreeVent PUR Basic Size#7/65mm	7	10,6	7,0	8,7	5,6	65	F19558070
FreeVent PUR Basic Size#8/70mm	8	12,5	8,3	10,5	6,1	70	F19558080
FreeVent PUR Basic Size#9/75mm	9	13,6	9,4	11,6	7,1	75	F19558090
FreeVent PUR Basic Size#10/85mm	10	14,6	10,3	12,5	8,0	85	F19558100
FreeVent PUR Basic Size#11/87mm	11	15,3	11,2	13,5	9,0	87	F19558110
FreeVent PUR Basic Size#12/90mm	12	16,5	12,4	14,5	10,2	90	F19558120



PUR Basic HME



22 mm HME connector



FreeVent PUR Basic HME Tracheostomy Set consists of:

- 1 outer tube
- 1 plain inner cannula
- 1 inner cannula with 22 mm HME connector

The PUR tracheostomy tubes provide a direct tracheal access for long term airway management. Thermosensitive PUR – provides sufficient rigidity for insertion but softens at body temperature which increases the comfort for the patient. This version has an inner cannula fitted with an HME-connector to simplify the use of an HME Heat and Moisture Exchanger). The unique conical shape also makes it much easier to insert.

Feature

Thermoplastic material
 Conical shape
 Wide range of models & sizes
 Low weight
 Excellent fit between inner and outer tube
 Low profile neck plate
 HME connector

Benefit

Soft and flexible at body temperature
 Easier insertion & comfort
 Fits all
 High comfort
 No leakage
 Discrete design
 Easier to use HME

Order information:

Product name	Size	Outer tube		Inner cannula		Length in mm	REF
		OD in mm	ID in mm	OD in mm	ID in mm		
FreeVent PUR Basic HME Size#7/65mm	7	10,6	7,0	8,7	5,6	65	F1511070
FreeVent PUR Basic HME Size#8/70mm	8	12,5	8,3	10,5	6,1	70	F1511080
FreeVent PUR Basic HME Size#9/75mm	9	13,6	9,4	11,6	7,1	75	F1511090
FreeVent PUR Basic HME Size#10/85mm	10	14,6	10,3	12,5	8,0	85	F1511100
FreeVent PUR Basic HME Size#11/87mm	11	15,3	11,2	13,5	9,0	87	F1511110
FreeVent PUR Basic HME Size#12/90mm	12	16,5	12,4	14,5	10,2	90	F1511120



Inner cannula with Speech Valve round

FreeVent PUR Speech Tracheostomy Set consists of:

- 1 outer tube with fenestration sieve
- 1 inner cannula with speech valve round window fenestration
- 1 inner cannula with 15 mm ISO connector with window fenestration

The PUR tracheostomy tubes provide a direct tracheal access for long term airway management. Thermosensitive PUR – provides sufficient rigidity for insertion but softens at body temperature which increases the comfort for the patient. This version has an inner cannula fitted with a Speech Valve to simplify the use of an HME (Heat and Moisture Exchanger). The unique conical shape also makes it much easier to insert.

Feature	Benefit
Thermoplastic material	Soft and flexible at body temperature
Conical shape	Easier insertion & comfort
Wide range of models & sizes	Fits all
Low weight	High comfort
Excellent fit between inner and outer tube	No leakage
Low profile neck plate	Discrete design
15mm ISO connector	Fits TrachPhone HME

Order information:

Product name	Size	Outer tube		Inner cannula		Length in mm	REF
		OD in mm	ID in mm	OD in mm	ID in mm		
FreeVent PUR Speech Size#7/65mm	7	10,6	7,0	8,7	5,6	65	F19511070
FreeVent PUR Speech Size#8/70mm	8	12,5	8,3	10,5	6,1	70	F19511080
FreeVent PUR Speech Size#9/75mm	9	13,6	9,4	11,6	7,1	75	F19511090
FreeVent PUR Speech Size#10/85mm	10	14,6	10,3	12,5	8,0	85	F19511100
FreeVent PUR Speech Size#11/87mm	11	15,3	11,2	13,5	9,0	87	F19511110
FreeVent PUR Speech Size#12/90mm	12	16,5	12,4	14,5	10,2	90	F19511120



Inner cannula with HME connector



FreeVent PUR Speech HME Tracheostomy Set consists of:

- 1 outer tube with fenestration sieve
- 1 inner cannula with speech round valve window fenestration
- 1 inner cannula with 22 mm HME connector with window fenestration

The PUR tracheostomy tubes provide a direct tracheal access for long term airway management. Thermosensitive PUR – provides sufficient rigidity for insertion but softens at body temperature which increases the comfort for the patient. The unique conical shape also makes it much easier to insert.

Feature	Benefit
Thermoplastic material	Soft and flexible at body temperature
Conical shape	Easier insertion & comfort
Wide range of models & sizes	Fits all
Low weight	High comfort
Excellent fit between inner and outer tube	No leakage
Low profile neck plate	Discrete design
HME connector	Easier to use HME

Order information:

Product name	Size	Outer tube		Inner cannula		Length in mm	REF
		OD in mm	ID in mm	OD in mm	ID in mm		
FreeVent PUR Speech HME Size#7/65mm	7	10,6	7,0	8,7	5,6	65	F1460070
FreeVent PUR Speech HME Size#8/70mm	8	12,5	8,3	10,5	6,1	70	F1460080
FreeVent PUR Speech HME Size#9/75mm	9	13,6	9,4	11,6	7,1	75	F1460090
FreeVent PUR Speech HME Size#10/85mm	10	14,6	10,3	12,5	8,0	85	F1460100
FreeVent PUR Speech HME Size#11/87mm	11	15,3	11,2	13,5	9,0	87	F1460110
FreeVent PUR Speech HME Size#12/90mm	12	16,5	12,4	14,5	10,2	90	F1460120

HME - Adapters

- Silicone adapter that fits inner cannula.

Order information:

Product name	Quantity	REF
HME adapter for FreeVent PUR	1	F19872

Order information:

Product name	Quantity	REF
HME adapter for FreeVent Silver	1	F19871



HME adapter for PUR and Silver

FreeVent Speech Valve

- Speech valve round that fits to PUR inner cannula AD.
- With or without O₂ connection.

Order information:

Product name	Quantity	REF
FreeVent speech valve	1	F19509



FreeVent speech valve

Order information:

Product name	Quantity	REF
FreeVent speech valve with O ₂ connector	1	F19563



FreeVent speech valve with O₂ connection

FreeVent Plug

- FreeVent plug is used to close the tube in order to wean the patient.

Order information:

Product name	Quantity	REF
FreeVent plug	1	F10682



FreeVent plug

FreeVent- spare parts and inner cannulas

All FreeVent spare parts and inner cannulas can be ordered separately i.e all connectors in both PUR or Silver material can be ordered. Contact your local dealer for more information.



Example: Extra inner cannula AD in PUR